

# Behind the ear

# Ahead of the curve



### Introducing the new Stride P

Your patients are going to love Stride P. The inspired design of this wireless BTE puts control at their fingertips with an integrated and intuitive push button and volume control.

### Stride at a glance

		Pro	800	700	600	500
Channels		20	20	16	10	6
Signature features		Glossary				
SpeechZone 2	In noisy situations, adjusts to an optimal microphone strategy to focus on speech from any direction	SpeechZone 2	SpeechZone			
Binaural spatial processing	Hearing instruments work together to determine the direction of speech	•	•			
SoundNav	Automatic program that identifies and smoothly transitions across environments, with a special focus on conversations	7 environments	6 environments	5 environments	2 environments	AutoMic
Sound Conductor	Dynamically balances the appropriate features providing speech understanding and comfort while maintaining natural sound quality	•	•	•	•	•
Log It All	Logs time spent across all seven SoundNav environments to support technology recommendations	•	•	•	•	•
MyMusic	Specialized fitter-controlled settings to enhance music enjoyment	Automatic binaural	Automatic binaural	•	•	•
Binaural Phone	Streams conversation from phone ear to opposite ear, no accessory required	•	•	•	•	
Automatic Adaptation Manager	Intelligently and gently transitions patients to the amplification required for long-term optimal speech understanding from first fit acceptance	•	•	•	•	•

### Features

Adaptive directional	Directional strategy that tracks and suppresses multiple background noise sources	Multiband	Multiband	Multiband	Multiband	•
Frequency compression	Compresses and shifts inaudible high-frequency sounds into a more audible lower frequency region	•	•	•	•	•
AntiShock	Eliminates discomfort from sudden impulse noises	•	•	•	•	•

### In all technology levels:

- Pinna Effect
  - Natural Sound Balance
  - Data logging
  - Feedback manager
  - Wind control
  - Tinnitus masker
  - Manual programs
  - Streaming programs
  - DuoLink
  - Easy-t
  - Easy-DAI
- IntelliVent technology for custom ear pieces
  - Plasma coating
  - IP67 for Stride P and Stride P Dura
  - Telecoil
  - WDRC processing with NAL-NL1, NAL-NL2 and DSL v5
  - Linear processing with DSL v5

Compatibility may be limited by product style

Accessories give you another tool to solve patient challenges. They are available at all technology levels.

unitron.com

© 2015 Unitron. All rights reserved 15-024 028-6176-02

Stride™  
Powered by North

unitron. Hearing matters



# Help every patient find their Stride

Stride™ leads the way in patient-centric design. This new family of behind-the-ear (BTE) and in-the-ear (ITE) styles offers the impressive aesthetics, comfortable fit and intuitive controls that patients find so appealing, with choices to satisfy every need. Since it's built on North, the platform that's powering a new direction in natural sound, Stride is packed with technology that truly allows your patients to focus on conversations in all of their listening environments. It's never been easier to help your patients find their Stride.



Stride BTEs

Stride P is a new wireless 13 battery BTE that pairs stylish design with dual controls, while Stride P Dura is designed with water resistance and durability in mind.




Stride ITEs

An all-new lineup of Stride ITEs are freshly designed inside and out, featuring a variety of custom products, intuitive controls, enhanced durability and optimal comfort.

# Step into conversations with confidence

Our new SoundNav automatic program offers precise classification of seven different environments, four of which are specifically focused on conversations. SoundNav works in harmony with other trailblazing technologies – Sound Conductor and SpeechZone™ 2 – to let patients seamlessly and automatically participate in conversations like never before, experience comfort with awareness in quiet and noisy environments, and enjoy the pure sound of music.

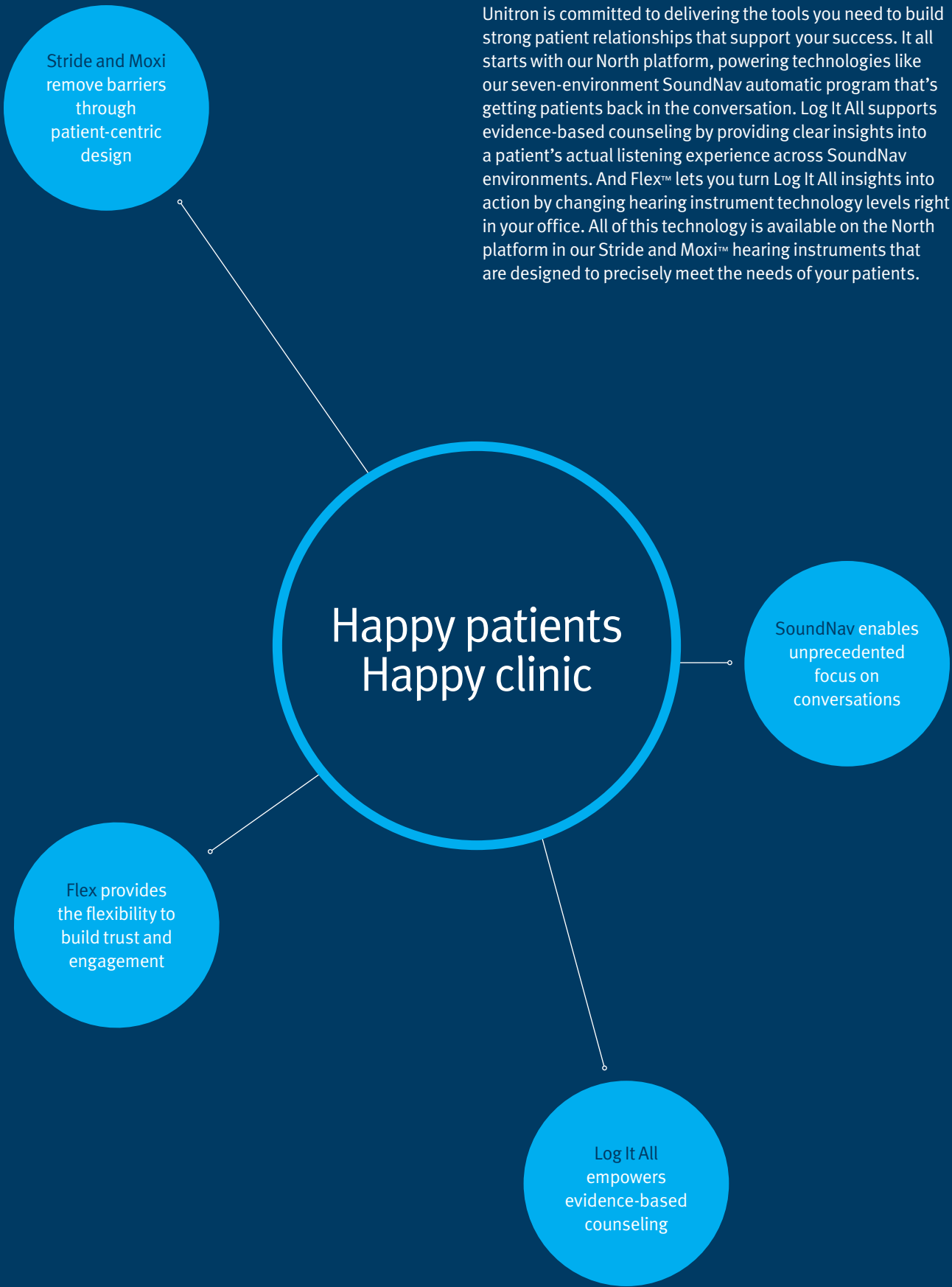

Unitron is committed to delivering the tools you need to build strong patient relationships that support your success. It all starts with our North platform, powering technologies like our seven-environment SoundNav automatic program that's getting patients back in the conversation. Log It All supports evidence-based counseling by providing clear insights into a patient's actual listening experience across SoundNav environments. And Flex™ lets you turn Log It All insights into action by changing hearing instrument technology levels right in your office. All of this technology is available on the North platform in our Stride and Moxi™ hearing instruments that are designed to precisely meet the needs of your patients.